


## Digital Waterbath

- Application: Biotechnology, Clinical, Environmental, Medical, Petroleum and Industrial etc.
- Concentric rings can be removed to accommodate containers of various sizes.
- Double walled chamber & Lid made of stainless steel
- The gap between the out and inner walls is filled with special grade glass wool insulation.
- Digital temperature controller (Optional)
- Temperature is controlled by capillary type thermostat.

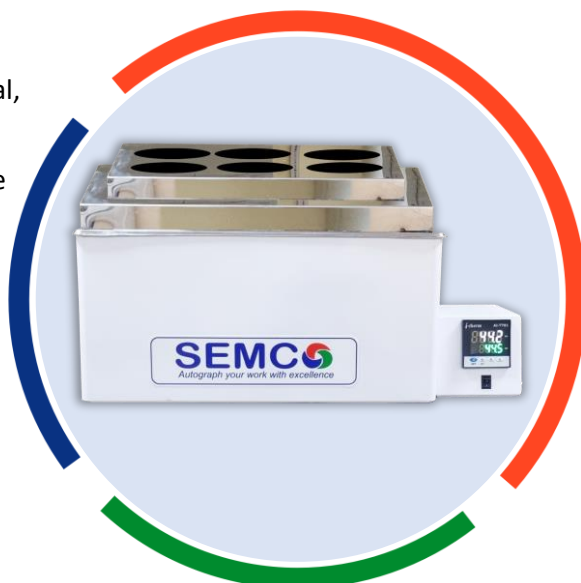

## Technical Specifications

Models		SSWB-10	SSWB-15	SSWB-20
Capacity & No. of Rings		10 Lit. / 6 Rings	15 Lit. / 9 rings	20 Lit. / 12 Rings
Temp.	Range	Ambient +5°C to 99°C		
	Accuracy	±1°C at 60°C		
	Uniformity	±2.0°C at 90°C		
	Sensor	PT100		
Heat-Up time		20 min. to 99°C		
Heating Power/Consumption		2kW	3kW	2kW
Display		Digital display (Optional)		
Controller		Capillary type thermostat / digital microprocessor based PID controller		
Lid		Stainless steel flat Lid		
Timer		Optional		
Safety Circuit		Over temp. and over current protection		
Drain Valve		Mounted		
Material		Internal: S.S., Exterior: Powder coated S.S., Insulation: Glass wool		
Power Supply		Single Phase, AC 220 V, 50Hz		

**NOTE: Technical specification & design are subject to change without prior notice due to continuous upgradation.**

### SEMCO SCIENTIFIC PRIVATE LIMITED

9C/91, Vrindavan Yojna, Raebareli Road, Near NKM School, Lucknow – 226029, U.P.  
 info@semcoscientific.in / semcoscientific@gmail.com / www.semcoscientific.in  
 +91-9918227744 / +91-6388224590 / 0522-3569875